VITALmaxx

Hornhaut-Entferner

06200 / 14016

DE

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ab Seite 3

FR

Mode d'emploi

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EN

Operating instructions

starting on page 10

NL

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vanaf pagina 24

INHALT

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CONTENTS

| | | | |
|-------------------------|----|----------------------|----|
| Items Supplied | 10 | Use | 14 |
| At A Glance | 10 | Cleaning and Storage | 15 |
| Symbols | 11 | Troubleshooting | 15 |
| Signal Words | 11 | Technical Data | 15 |
| Intended Use | 11 | Disposal | 16 |
| Safety Notices | 12 | Customer Service | 16 |
| Power Supply | 13 | | |
| Changing the Attachment | 14 | | |

ITEMS SUPPLIED

- 1 x basic device
- 1 x protective cap
- 2 x attachment (1 x coarse and 1 x fine)
- 1 x ring
- 1 x filter
- 1 x USB cable (A connector to C connector)
- 1 x cleaning brush
- 1 x operating instructions


Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department.

Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**

Four batteries (1.5 V, AAA type) are required for battery operation. These are not supplied.

AT A GLANCE

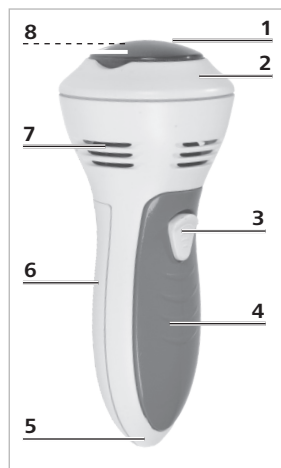- 1 Protective cap
- 2 Ring
- 3 On / Off switch
- 4 Handle
- 5 USB C input
- 6 Battery compartment
- 7 Ventilation slots
- 8 Attachment (under protective cap)

Dear Customer,

We are delighted that you have chosen the **VITALmaxx callus remover**.

Before using the device for the first time, please read the operating instructions carefully and keep them for future reference. If the device is given to someone else, it should always be accompanied by these operating instructions. They are an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

If you have any questions about the device and about spare parts/accessories, please contact the customer service department via our website: **www.ds-group.de/kundenservice**

SYMBOLS

Danger symbols: These symbols indicate possible dangers. Read the associated safety notices carefully and follow them.

Supplementary information

Read operating instructions before use!

Direct current

SIGNAL WORDS

Signal words indicate hazards if the associated notices are not followed.

DANGER – high risk, results in serious injury or death

WARNING – medium risk, may result in serious injury or death

CAUTION – low risk, may result in minor or moderate injury

NOTICE – may result in risk of damage to material.

INTENDED USE

- This device is intended to be used to remove calluses as well as dry skin from your feet.
- Only use the device on your feet; do not use it on other parts of your body.
- The device is designed for personal use only and is not intended for commercial applications.
- The device should only be used as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

SAFETY NOTICES

WARNING – Danger of Injury!

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- The device is **not** a toy!
- The device must always be switched off and disconnected from the mains power when it is left unattended and before assembling, disassembling, changing battery or cleaning the device.
- Keep the packaging material away from children and animals. There is a danger of suffocation.
- Make sure that the USB cable is always kept out of the reach of small children and animals. There is a danger of strangulation.
- Batteries can be fatal if swallowed. Keep batteries and the device out of the reach of children and animals. If a battery is swallowed, immediate medical attention must be sought. There is a danger of burns from battery acid.
- If a battery leaks, avoid physical contact with the battery acid. Wear protective gloves if necessary. If you touch battery acid, rinse the affected areas with abundant fresh water straight away and seek medical assistance immediately. There is a danger of burns from battery acid.
- People who have a pacemaker or belong to another risk group should consult a doctor before using the device.
- The device is not suitable for people who suffer from haemophilia, diabetes, or an immunodeficiency disorder.
- Do not use the device on injured, infected, loose, or (sun-)burnt skin or in the

direct vicinity of acne, birthmarks, warts or varicose veins.

- If you have skin problems on your feet, consult a doctor before using the device.
- For reasons of hygiene, an attachment should only be used by one person.
- Do not use the device if it is damaged. Before each use, check whether the device is in good condition. If there is any damage or wear or sharp edges, injuries could occur.
- Make sure that the rotating parts of the device only come into contact with the calluses to be treated. Never touch the attachment with your fingers during operation!
- Keep clothing, jewellery, hair, hands and other body parts away from the device so that nothing can get caught in rotating device parts. There is a danger of injury due to rotating device parts.
- Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.
- Consult a doctor if skin irritations occur.
- Make sure that the connected USB cable does not present a trip hazard.

DANGER – Danger of Electric Shock!

- No liquid must be allowed to enter the device. Do not immerse in water! Never use while bathing or showering!
- Never immerse the device and the USB connector in water or other liquids and ensure that they cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull it out of the water while it is still connected to the mains power!
- Never touch the device and the USB connector with wet hands.

NOTICE – Risk of Damage to Material and Property!

- Only use suitable USB mains adapters whose technical data is the same as that of the device.
- Make sure that the USB cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces. Do not wrap the USB cable around the device (danger of the cable breaking).
- The device should only be used if it has been fully and correctly assembled! **The filter must be inserted and must not be damaged!** If the dirt were to get inside the device, this could damage the motor.
- Do not insert any objects into the openings of the device.
- Protect the device from heat, flames, persistent direct sunlight, sub-zero temperatures, persistent moisture and impacts.
- Switch off the device and put on the protective cap when you are not using it.
- In the event of damage / faults, switch the device off immediately.
- Do not make any modifications to the device. Do not try to repair it yourself or to replace any individual parts. In the event of damage, contact the customer service department.
- Use only the battery type specified in the technical data. When replacing batteries, ensure the correct polarity (+/-).
- Always replace all batteries together. Never use batteries of different types or brands, or batteries with different capacities.
- Take the batteries out of the device when they are dead or if you will not be using the device for a prolonged period. This will prevent any damage that might occur if a battery leaks.
- Remove a dead battery from the device immediately.

POWER SUPPLY

Please Note!

- No batteries may be inserted in the device if it is to be operated via the USB cable.

The device can either be operated with batteries or be connected to the mains power via the USB cable.

Batteries

1. Pull up the battery compartment cover by the tab labelled **OPEN** and take it off.
2. Remove any used batteries if necessary and dispose of them in an environmentally friendly way.
3. Insert four new batteries. Pay attention to the correct polarity (+/-).
4. Insert the locking lug of the battery compartment cover into the opening which is provided for it and press the battery compartment cover shut. It must engage audibly.

Replace the batteries when the power of the device starts to wane.

USB Cable

1. Connect the end with the USB C connector to the USB C input (5) on the device.
2. Connect the other end with the USB A connector to a plug socket via a suitable USB mains adapter (not supplied with the device).

CHANGING THE ATTACHMENT

Please Note!

- A damaged / worn attachment can cause injuries to the skin and must be replaced.
1. Grip the protective cap **(1)** at the level of the two recesses on the ring **(2)** and pull it off.
 2. Grip the attachment **(8)** at the level of the two recesses on the ring and pull it off the device.
 3. Place another or new attachment on the mount and press it on firmly. It must engage audibly.

USE

Please Note!

- Clean the attachments before you first use them as they may have possible production residues on them. Follow the instructions in the "Cleaning and Storage" chapter to do this.
- The device should only be used with a correctly applied attachment!
- Only use the device on uninjured calluses on your feet.
- Do not use the device on tender, soft skin.
- Do not press the attachment against the skin with force. It could be injured.
- From time to time, check the area you are treating to avoid penetrating your skin too deeply and causing injuries. Do not hold the attachment on one spot for longer than 2 to 3 seconds.
- If you experience an unpleasant sensation while you are using the device, stop applying it immediately!

Removal of Calluses

For dry application, make sure that there are no cream or care oil residues on your feet.

1. Connect the device to the mains power via the USB cable or insert batteries (see the "Power Supply" chapter).
2. Grip the protective cap **(1)** at the level of the two recesses on the ring **(2)** and pull it off.
3. Place the required attachment **(8)** on the device (see the "Changing the Attachment" chapter).
4. Slide the On / Off switch **(3)** to the **ON** position to switch on the device. The attachment starts to rotate and the suction function is activated.
5. Place the rotating attachment against the callus and move the device back and forth and to the sides. Move the device slowly and without pressure.

If the device is pressed too firmly against the skin, the attachment will automatically come to a stop to prevent skin injuries.

6. Switch off the device after use by sliding the On / Off switch to the **OFF** position.

After Treatment

1. Apply cream to the feet so that the skin becomes soft and tender.
2. Clean the device after each use (see the "Cleaning and Storage" chapter).

CLEANING AND STORAGE

Please Note!

- Do not use any caustic or abrasive cleaning agents to clean the device. These may damage the surface.

1. Remove the attachment (8) from the device (see the “Changing the Attachment” chapter) and brush off any flakes of skin with the cleaning brush that is also supplied. Then rinse off the attachment under running water.
2. Turn the ring (2) towards the  symbol and take it off.
3. Remove the filter and tap it out. Clean it with the cleaning brush if necessary.
4. Place the filter on the device.
5. Place the ring on the device and turn it towards the  symbol to fix it in place.
6. Place the attachment on the mount in the device. **Please Note:** The attachment must be completely dry!
7. Place on the protective cap (1) and press it on firmly so that its locking lugs hook into the suction openings in the ring.
8. Wipe down the device with a damp cloth and then dry it with a dry cloth.
9. Store the device and accessories in a dry place which is not accessible to children and animals.

TECHNICAL DATA

| | |
|-------------------------------|--------------------------------------|
| Article number: | 06200; 14016 |
| Model: | BG-F-900 |
| Voltage supply: | 6 V DC (4 x 1.5 V battery, AAA type) |
| ID of operating instructions: | Z 06200_14016 M DS V5 0824 as |

TROUBLESHOOTING

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact the customer service department.

Do not attempt to repair an electrical device yourself!

The device does not work.

Battery mode:

The batteries are too weak or incorrectly inserted.

- ▶ Insert new batteries and check the polarity.

Mains operation:

The plug connections are not fitted correctly.

- ▶ Correct the fit.

The plug socket is defective.

- ▶ Try another plug socket.

The mains adapter is defective or unsuitable.

- ▶ Try another mains adapter.

The calluses are not being properly removed.

The attachment (8) is worn.

- ▶ Replace the attachment.

The device does not suck properly.

The filter is clogged up.

- ▶ Clean the filter. The filter must be replaced if necessary.

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

(Rechargeable) batteries must not be disposed of along with the household waste. Consumers are legally obliged to dispose of (rechargeable) batteries separately.

(Rechargeable) batteries can be handed in free of charge at a collection point of the local authority/district or in retailers so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants.

Only dispose of (rechargeable) batteries when they are discharged. If possible, use rechargeable batteries rather than single-use batteries.

In the case of (rechargeable) batteries containing lithium, tape over the poles before you dispose of them to prevent a short circuit. A short circuit may cause a fire or explosion.

The adjacent symbol (crossed-out bin underlined) means that used devices should not be disposed of in the domestic waste, but should instead be taken to special collection and return systems.

Owners of used devices must remove from the used device any used batteries or rechargeable batteries which are not encapsulated by the used device and can be removed in a non-destructive way and dispose of them separately (see also the section on battery disposal).

Owners of used devices from private households can take them free of charge to the collection points of the public waste disposal authorities or the disposal points which are set up by manufacturers and distributors under the Electrical Equipment Act so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants. Businesses which market electrical and electronic devices are also obliged to take back used devices.

Waste prevention makes an even more valuable contribution to protecting the environment. This means if possible, as well as continuing to use or repair a device yourself, passing it on to another user also represents an ecologically beneficial alternative to disposal.

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Stormarnring 14, 22145 Stapelfeld, Germany

CUSTOMER SERVICE

Am Heisterbusch 1, 19258 Gallin, Germany

☎ +49 38851 314650* | kundenservice@dsp.de

* Calls to German landlines are subject to your provider's charges.

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All rights reserved.
DS Produkte GmbH