A woman's back is shown from the side, with her hair tied back and a purple flower in it. Overlaid on her back are four red heat therapy symbols, each consisting of a central red circle with three concentric red arcs radiating outwards. The background features a light gray geometric pattern of overlapping squares.

Backnetix[®]


WRAP

WELLNESS




International patents pending. All imitations will be prosecuted.
BACKNETIX[®] WRAP is a registered EU/CTM trade mark.

Attaching system for neck
 Sistema de sujeción para cuello
 Système de fixation pour le cou
 Befestigungssystem im Halsbereich
 Sistema di fissaggio per il collo
 Sistema de fixação para pescoço
 Bevestigingssysteem voor de nek

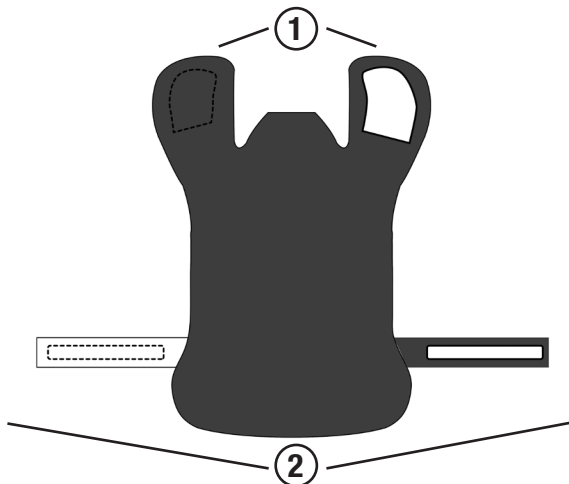Attaching system for waist
 Sistema de sujeción para cintura
 Système de fixation pour la ceinture
 Befestigungssystem im Hüftbereich
 Sistema di fissaggio per la vita
 Sistema de fixação para cintura
 Bevestigingssysteem voor de taille

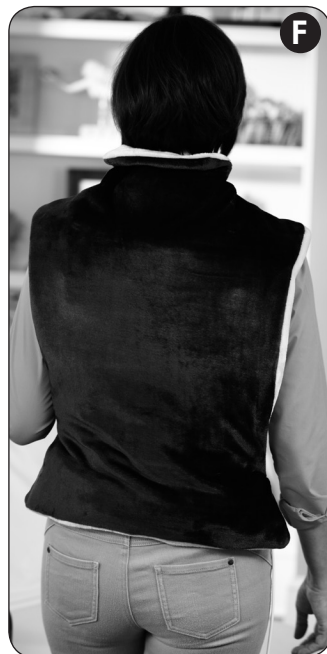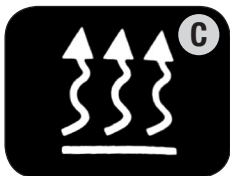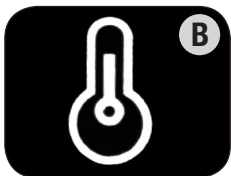Thank you for purchasing BACKNETIX® WRAP, the neck and back wrap that applies heat and massage to help relieve tension.

Before using the product for the first time, please read all these instructions so you use the product correctly.

1. WARNINGS

Read the instructions carefully before using the device for the first time. Store them for future reference.

It is absolutely essential for the user not to fall asleep while wearing the electric wrap.

If the wrap is worn for long time periods, the temperature selector should be set to the continuous position.

Do not use the electric wrap in a folded position! Do not scrunch the wrap up!

Do not stick pins into the electric wrap!

Do not use the electric wrap on animals!

Electric Shock Hazard!

Connect and use the device only if the mains supply complies with the information on the device's rating plate.

Check the electric wrap regularly to see if it shows signs of wear or damage. Do not connect the device to the mains if the power cable, temperature selector, plug or the wrap itself shows any visible signs of damage.

If the device is broken, unplug it immediately.

To prevent hazardous situations, any repairs or work on the device should be done by a qualified technician from the brand's official technical support service.

Check the device's power cable regularly:

Make sure it never makes contact with hot objects or parts.

Do not expose it to sharp edges or corners.

Use the electric wrap **only in dry rooms** (i.e. do not use it in the bathroom or similar places).

Do not use the electric wrap **if it is wet or damp!** Ensure the device is completely dry before using it.

The wrap creates low-intensity electrical and magnetic fields that can have a negative effect on pacemaker functions. We therefore recommend contacting your doctor or the manufacturer of the electric wrap before using it.

Danger of Overheating!

Do not use the electric wrap **without the supervision** of someone who knows how to use the device!

Do not cover the temperature selector.

To ensure the temperature selector works properly, it should never be placed on top of or under the electric wrap. Nor should it be covered by any other object.

Never wrap the power cord around the temperature selector!

Continuous use of the electric wrap at high temperatures may cause burns to the skin.

Danger of Injury!

The electric wrap should not be used by the disabled, people who need assistance or those who are insensitive to heat. It should not be used by young children.

Do not allow persons with restricted physical sensory perception or mental abilities or with lack of experience and knowledge to operate the device unless they are supervised or have been instructed on the use of the appliance by somebody who is responsible for their safety. Keep out of reach of children. Supervise children to prevent them from playing with the device.

When storing the appliance, allow it to cool down before folding;

Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again;

This appliance is not intended for medical use in hospitals;

Children under the age of three are not to use this appliance due to their inability to react to overheating.

The appliance is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or unless the child has been adequately instructed on how to operate the controls safely.

The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

Do not stick pins into the electric wrap!

Do not use the electric wrap in a folded position.

Not suitable for children under 3 years of age.

Do not wash.

DOUBLE INSULATION

This appliance is double insulated; therefore, no earth wire is required.

Always check that the mains voltage corresponds to the voltage on the product rating plate.

2. THE SET INCLUDES:

- 1 BACKNETIX® WRAP
- 1 Instruction manual
- 1 Travel bag

3. INSTRUCTIONS FOR USE

The BACKNETIX® WRAP is specifically designed to fit around the neck, shoulders, and upper and lower back.

1. Remove the BACKNETIX® WRAP from its protective plastic packaging and unfold it.
2. Place the BACKNETIX® WRAP on your shoulders and back with the blue side up. Make sure it fits around the neck, as shown in *Picture F*.
3. Adjust the attaching systems at the neck and waist to ensure the BACKNETIX® WRAP is in full contact with the body.

Picture 1: Attaching system for neck.

Picture 2: Attaching system for waist.

4. The BACKNETIX® WRAP has two heat levels and four massage levels that can be used individually or in combination for customized therapy (*Figure A*).

5. Plug the device into the mains. The red pilot light on the lower part of the temperature selector will come on when the device is plugged into the mains.

To activate the BACKNETIX® WRAP, press the temperature button (*Picture B*) and/or the massage button (*Picture C*) to activate their respective functions.

6. Press the temperature button (*Picture B*) to select the desired heat level. The green pilot light indicates low temperature and the red pilot light indicates high temperature.

7. Press the massage button (*Picture C*) to select the vibration intensity and vibration mode:

1. Low vibration – constant – green pilot light
2. Low vibration – intermittent – flashing green pilot light
3. High vibration – constant – red pilot light
4. High vibration – intermittent – flashing red pilot light

Note: The heat function and massage function can be operated individually or at the same time.

8. To turn off the wrap, press the temperature and/or massage buttons until the device shuts off.

NOTE: The pilot light on the lower part of the temperature selector will start flashing if there is any kind of product error. If this happens, unplug the product from the mains and check the power cable and plug.

Burn Hazard

When wearing the electric wrap, it should not be too tight, but should be worn loose against the body. Never wear the device without the outer case!

4. AUTOMATIC SHUTOFF FUNCTION

The BACKNETIX® WRAP will shut off automatically after two hours of continuous use of the temperature function and after 30 minutes of continuous use of the massage function.

5. TECHNICAL SPECIFICATIONS

Voltage: 220-240V

Frequency: 50 Hz

Power: 100 W

Material of outer case: 100% polyester

Material of inner case: 100% polyester

Material of handset: ABS

6. CLEANING

Electric shock hazard!

Before cleaning the wrap, unplug the power cable from the mains!

IMPORTANT: Only the outer case is washable. Do not wash the inner case of the BACKNETIX® WRAP.

(*Figure D*) Outer case

(*Figure E*) Inner case

The outer case cannot be dry cleaned.

Do not use bleach on the outer case.

Do not wring dry the outer case.

Do not get the inner case (connected to the power cable) damp or wet.

Unplug the BACKNETIX® WRAP from the mains and remove the outer case.

Remove the outer case by undoing the zip on the side.

The outer case is machine washable in warm water on gentle cycle. Tumble dry on medium heat. Remove at once when tumbling stops.

The outer case can also be hand washed in warm water. Air dry.

When the outer case is completely dry, insert the inner case inside to use the wrap again.

When not in use, store method/way/place:

- Keep the BACKNETIX WRAP clean and dry.
- Store the device in a dry place when it is not in use.
- Carefully fold the product to store it.
- Allow the device to cool before you fold or roll it up.
- Please do not use any chemical moth protection.
- Do not place any heavy objects on the device during storage, as this may damage the electrical wiring.

CAUTION

DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations. Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

QUALITY WARRANTY

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993, both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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ISL GmbH
Best Direct B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.