# promed ACP

Elektro-Akupunkturgerät  
Electric Acupuncture Device  
Appareil d'acupuncture électronique  
Apparecchio per elettroagopuntura  
Dispositivo eléctrico de acupuntura  
Elektrisch acupunctuurtoestel  
Aparat do elektroakupunktury  
Прибор для электроakupunktury  
Elektrikli Akupunktur Cihazı

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

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DESCRIPTION OF THE DEVICE	p.	18	MANUFACTURER	p.	144
SCOPE OF DELIVERY / SPARE PARTS	p.	19	TYPE PLATE	p.	144
SAFETY INFORMATION	p.	20	CERTIFICATION NUMBER	p.	144
THINGS TO KNOW ABOUT ACUPUNCTURE	p.	21	NORMALIZED SYMBOLS	p.	128
OPERATION OF THE DEVICE	p.	22	CE DIRECTIVES	p.	130
FREQUENTLY ASKED QUESTIONS AND ANSWERS	p.	25	TECHNICAL SPECIFICATIONS	p.	132
TROUBLESHOOTING	p.	27	TRANSPORT/STORAGE/OPERATING CONDITIONS	p.	134
CARE, MAINTENANCE AND STORAGE	p.	28	APPLICATION EXAMPLES	p.	136
DISPOSAL	p.	28	GUARANTEE CARD	p.	145
GUARANTEE	p.	29			

DE

EN

FR

IT

ES

NL

PL

RU

TR


AR

Please read the user manual thoroughly before using the device for the first time.

Thank you for purchasing the **Promed ACP**. You have purchased a high-quality product designed for personal hygiene and wellbeing purposes. **Promed** is a leading company with decades of experience in personal hygiene, wellness and health areas. The **Promed ACP** device has been designed and manufactured in compliance with Medical Devices Directive 93/42/EEC in order to guarantee safe use and may be used after reading these operating instructions. We, as the manufacturer, cannot be made liable in any way for injury or damages to people or objects that arise from failure to comply with this user manual. We wish you a lot of enjoyment with your new **Promed ACP**. We would like to familiarize you with the therapie unit in the following sections. Please read the user manual thoroughly before using the device for the first time.

You can receive professional advice wherever **Promed** products are sold or get in touch with us if you have any questions. We can provide you with the name of the representative responsible for you.

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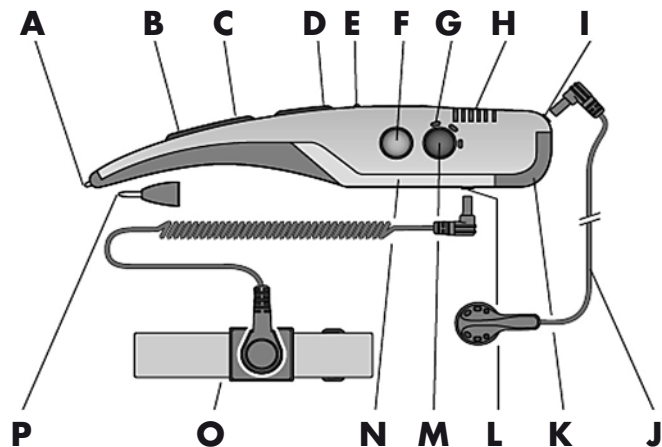Lindenweg 11, D-82490 Farchant, Ph: +49 (0) 8821/9621-0,

Fax: +49 (0) 8821/9621-21, E-Mail: info@promed.de

Further informations can be found on our homepage:

**www.promed.de**

- A.** Acupuncture tip
- B.** Display "Stimulation"
- C.** Stimulation button
- D.** Intensity button (-/+)
- E.** Battery display
- F.** Power on / off button
- G.** Time display (manual, 30 s, 60 s)
- H.** Intensity display / Resistance display
- I.** Earphone connection
- J.** Earphone
- K.** Battery compartment lid
- L.** External electrode connection
- M.** Time selector button
- N.** Contact electrode
- O.** External electrode with cable
- P.** Acupuncture tip for ear acupuncture

The following is included in the scope of delivery of your

**Promed ACP** kit:

- 1x electric acupuncture device
- 1x acupuncture tip
- 2x CR 2032 button cells, already inserted
- 1x earphone
- 1x external electrode with cable
- 2x companion books (1x German, 1x English)
- 1x operating manual

We reserve the right to apply changes for the purposes of technical improvement without prior notice.

The following spare parts can be purchased from **Promed**:

- Acupuncture tip (article no.: 440011)
- Earphone (article no.: 440012)
- External electrode with cable (article no.: 440013)

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

**Please consider before commissioning:**

- Only use the device for the intended purpose as specified by the operating instructions.
- Only for use in indoor spaces.
- Children must be supervised to ensure that they do not play with the device.
- Keep the device away from children and pets.
- When using a device in the presence of children, attentive supervision is required.
- Never position or use the device in a wet or moist environment.
- Do not use the device under water, e.g. in the shower.
- Never bring the device into direct contact with fire, gas or oxygen, nor with hot objects, such as hobs.
- Take every possible precaution to ensure that the device does not fall down or is otherwise damaged.
- If problems with the device occur, please have it repaired.
- Do not lubricate or wash the device.

**Danger!**

- Do not bring the device into direct contact with water or other liquids. Should any liquid enter the device, immediately remove the batteries and refrain from further use. Contact your dealer or our service department.
- The device must not be used outdoors.
- Never touch the device with wet hands.

- Do not store the device near a sink or a bathtub as it may fall or be pulled into the sink or the bathtub.
- Never position the cables around your neck, to prevent any risk of strangulation.
- Do not put any small parts into your mouth, to prevent inhaling or swallowing them.

**Warning!**

- The device is not intended for use with pets or children.
- Do not leave the device without supervision when there are children or persons not used to handling the device around.
- Do not use any spare parts not recommended by the manufacturer. Only use the external electrode supplied with the device. As the manufacturer, we do not assume any responsibility if you use another electrode and this does not work, the device is damaged or if this causes an uncomfortable sensation and/or injury during application.
- Ensure that children do not play with the device.
- This device is not intended for use by persons (including children) with limited physical, sensory or mental capabilities or lacking the required experience and/or knowledge, except if they are supervised by a person responsible for their safety and have been instructed by said person on how to use the device.
- The device must not be used outdoors.
- Do not allow children to play with the packing material as there is a risk of suffocation.

- Do not insert any foreign objects into any of the openings. Avoid any contact between the device and pointed or sharp objects. Using the device for anything but its intended purpose will void the warranty.
- In case of diabetes or other illnesses, please consult your general practitioner before using the device.
- Pregnant women should take the necessary precautions and consider their individual resilience; consult your doctor if necessary.
- Do not treat any body parts showing signs of swelling, burns, inflammation, skin rashes, wounds or abnormal sensitivity.
- The treatment should be comfortable. If you feel pain or experience the application as uncomfortable, discontinue use and consult your doctor.
- Therapy with **Promed ACP** is not a replacement for medical diagnoses or treatments. Consult your doctor in any case of pain or illness prior to using the device.
- Do not use the device if you have a pacemaker or any other implants.
- Electrical stimulation must not be applied while operating machinery or while driving.
- If any irregularities of the heart or cardiac diseases are known (e.g. arrhythmia or myocardial damage), the device should not be used.
- If electrodes are used, they must not be applied to open skin parts. Remove the electrodes prior to bathing or swimming.
- The electrical acupuncture tip and the additional tip are made of gold-plated brass. In rare cases, allergic skin reactions may occur.
- If skin irritation occurs, discontinue treatment and consult your doctor.
- The device must not be applied above or through the head (transcerebrally), directly to the eyes, above the mouth, on the front of the neck, near the larynx or across the heart (neither in the chest or back region nor through the upper body).
- Acupuncture therapy should only be applied to the points listed in the companion book. For additional information, ask your doctor or natural health professional.

### THINGS TO KNOW ABOUT ACUPUNCTURE

Acupuncture is a traditional Chinese form of therapy in which specifically defined body parts, the so-called acupuncture points, are stimulated by needles. They are located along body lines called meridians. These acupuncture points can be used to balance energetic disorders in the body, stimulate specific organs or body regions and foster their function and vitality. This kind of stimulation can also be obtained by applying pressure (acupressure) or harmless electric impulses, like with the **Promed ACP** electric acupuncture device, to the corresponding acupuncture points. Acupuncture is used to treat pain syndromes and functional disorders, among other things.

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

### Initial commissioning

1. Remove the battery compartment lid (**K**) on the underside of the device by sliding it backwards.
2. Pull out the contact protection film from between the batteries.
3. Close the battery compartment.

### Testing / replacing the batteries

1. Once the battery power has been used up, the red LED of the battery display (**E**) starts flashing. To ensure proper operation, please replace the batteries in such a case.
2. To do so, remove the battery compartment lid (**K**) on the underside of the device by sliding it backwards.
3. Remove the two batteries and replace them with two new batteries of the same type (2 x 3 VDC lithium CR 2032). When inserting the batteries, ensure that the + sides point towards the underside of the device (the side with the rating plate).
4. Close the battery compartment.

#### Caution:

1. Remove the batteries when the device is not used for an extended period of time.
2. Do not mix old and new batteries or different battery types.
3. **Warning:** Should batteries leak and come into contact with skin or eyes, rinse these areas thoroughly with plenty of water.
4. Batteries should be used by adults only. Keep batteries out of reach of children.
5. Only batteries of the same or a similar type are recommended.

6. Remove empty batteries from the device.
7. Dispose of the batteries in accordance with the instructions provided by the battery manufacturer.

### Fields of application

The **Promed ACP** device can be applied in three fields:

- **Detection mode:** This is used to detect acupuncture points.
- **Manual therapy:** No time limit for therapy.
- **Therapy with automatic time settings:** 30 seconds or 60 seconds of therapy time for the acupuncture point.

Pick up the device as shown in the image. Hold the **Promed ACP** device in such a way that you are touching the contact electrode (**N**). A closed circuit is required for application. The contact electrode (**N**) forms one pole (counter-electrode) and the acupuncture tip (**A**) forms the other pole (active electrode). When applying treatment to another person, it is equally important to keep the circuit closed, i.e. by touching the person to be treated using the free hand or by using the external electrode (**O**).

### Detection mode

Switch on the device by pressing the Power On / Off button (**F**). An acoustic signal and a lit LED (**G**) above the time interval button (**M**) indicate that the device is switched on. The device is now in detection mode.

Apply the acupuncture tip (**A**) to the skin using slight pressure – approximately in the area where the acupuncture point should be located (see companion book). When you start the detection, the acoustic signal will become higher and one or more LEDs will

light up on the display **(H)** the nearer you get to an acupuncture point. Use the tip to look for the point until the acoustic signal is highest and a maximum of LEDs light up.

### PLEASE NOTE

Occasionally, acupuncture points may be located very close to each other. To find the proper point for the therapy you desire, use the application examples in the back of the instruction manual or in the companion book for orientation, or use additional literature on acupuncture.

### Manual therapy

Once you have found the acupuncture point, press the stimulation button **(C)** using the index or middle finger and hold. The acoustic signal will stop and the LED **(B)** will start to flash; the first LED on the display **(H)** will always flash simultaneously. You are now treating the acupuncture point with the lowest amperage. If you want to increase the intensity, we recommend keeping the acupuncture tip on the point and changing the intensity settings **(D)**. As soon as you let go of the stimulation button **(C)**, stimulation will stop and the device will return to the detection mode. This means you determine yourself how long the electrical impulses are sent.

### Therapy with automatic time settings

An internal clock can be used to set the **Promed ACP** device to 30 seconds or 60 seconds of stimulation time via the time selector button **(M)**. This new time setting will be maintained until it is reset by pressing the time selector button **(M)** again or switching off the device. In the basic settings, when switching on the device, the manual mode is always enabled initially. Prior to the detection of the acupuncture point, press the time selector button **(M)** to select the desired time: press once for 30 seconds and again for 60 seconds. The LEDs **(G)** indicate which setting has been selected. Detect the acupuncture point and briefly press the stimulation button **(C)** once. The LED **(B)** and the first LED on the display **(H)** will start to flash. Since you do not need to keep the stimulation button **(C)** pressed, you can comfortably apply slight pressure to the acupuncture point with the device. Once the set time has elapsed, the device will return to the detection mode. As the time settings are maintained, you can stimulate the point for the same pre-set time by again pressing the stimulation button **(C)** briefly.

### Changing the intensity

When switching on the **Promed ACP** device, it is always set to the lowest intensity level. For certain acupuncture applications, however, higher intensities may be required. You can change the application intensity prior to or during stimulation. It is recommended, however, to adjust the intensity only during stimulation, as different regions of the body may display different sensitivities and some people may be more sensitive to electrical current than others. With the device switched on, press the

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

intensity button **(D)** in the place labelled with a plus sign (+) several times, until you have reached the desired intensity. On the display, you can see that the visual indicator **(H)** changes to the next level each time the button is pressed. Once the last LED lights up, the highest intensity has been reached. Press the intensity button in the place labelled with a minus sign (-) in order to reduce the intensity step by step all the way down to the lowest intensity.

 **CAUTION**

The feeling of intensity may vary between different regions of the body. Therefore, it is advisable to always lower the intensity to the lowest level prior to changing the acupuncture point!

 **PLEASE NOTE**

As conductivity may vary between different regions of the body, it may be necessary to adjust the amount of pressure applied by the acupuncture tip to the selected acupuncture point.

 **PLEASE NOTE**

Do not forget to switch off the device after use. In case you have used the earphone and have put it away together with the device, only the lit LED will indicate that the device is still on. Battery use will increase if you keep the **Promed ACP** device switched on while not using it.

### Use of device accessories

The **Promed ACP** device is supplied with an acupuncture tip, an earphone and an external electrode with cable.

#### Acupuncture tip

The additional acupuncture tip **(P)** features a finer end piece with a very narrow tip, permitting even more efficient stimulation of the acupuncture point. This may be useful when detecting, finding and stimulating points located very closely to each other, e.g. in the ear. Simply attach the tip to the acupuncture tip **(A)** of the device. Ensure a tight fit and contact with the device.

#### Earphone

The earphone **(J)** is plugged into the socket **(I)** for connection with the **Promed ACP** device. Now the acoustic signal will only be heard on the earphone.

#### External electrode with cable

Plug the external electrode **(O)** into the connection socket **(L)**. For individual therapy, tie the electrode around your own wrist or arm; when treating another person, tie it around the person's wrist or arm instead. Now you can hold the **Promed ACP** device as you wish without the need to ensure contact with the contact electrode or - when treating a partner - additionally having to touch the other person. The closed circuit is ensured by attaching the external electrode. This provides you with more freedom during treatment.

 **PLEASE NOTE**

If necessary, you can improve contact with the external electrode by slightly moistening the electrode prior to applying it to the skin.

Only use the external electrode supplied with the device. As the manufacturer, we do not assume any responsibility if you use another electrode and this does not work, the device is damaged or if this causes an uncomfortable sensation and/or injury during application.

- **Why does my device sometimes not report any acupuncture points?**

**1. Individual treatment**

Have you established contact with the contact electrode?	The circuit must be closed at all times. Hold the <b>Promed ACP</b> device so that at least one finger touches the contact electrode <b>(N)</b> .
--	---

**2. Application to partner**

Have you closed the circuit by touching your partner?	Ensure here also that the circuit is closed at all times.
For treatment with external electrode application: Is the electrode applied correctly, i.e. fitting tightly and securely on the grease-free skin?	
Is the connection cable firmly plugged into the connection socket?	

**3. Treatment using the acupuncture tip**

Is the attached tip in contact with the device?	If necessary, remove it briefly, then reattach it and turn it slightly in order to re-establish contact.
---	--

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

- **How do I turn off the sound?**

Plug the earphone **(J)** into the connection socket **(I)**. Now you can no longer hear the sound in the application environment, but only directly over the earphone. In addition, the visual signals still provide further indication of the acupuncture points.

- **Can the current harm me?**

No, the amperage of the electrical impulses is so low that it will remain completely harmless even during longer application. However, the intensity should be adjusted depending on the sensitivity of the application point. For example, the skin in the facial region is more sensitive than the skin on the hand. For this reason, always start at the lowest intensity level for the stimulation of each new point and slowly increase it depending on individual preference. The application should never feel uncomfortable.

- **Which intensity should I choose?**

Please regulate the intensity in a way that it feels comfortable. It should never feel painful; however, you should be able to feel the impulses. The intensity will also depend on the location of acupuncture.

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

Issue	Potential cause	Solution
LEDs do not light up	No battery contact	Insert new batteries.
		Ensure that the batteries have been inserted correctly. Inspect the contacts for the following: <ul style="list-style-type: none"> <li>• All contacts are in the right place</li> <li>• No contacts have been damaged</li> </ul>
Stimulation is weak or cannot be felt	The acupuncture tip is not positioned correctly	Reposition acupuncture tip
	Intensity too low	Increase intensity level
	The batteries are empty	Replace with new batteries
Stimulation is uncomfortable	Intensity too high	Decrease intensity
	The device has not been operated as specified by the manual	Read the manual before use
Stimulation ineffective	Incorrect placement of acupuncture tip	Reposition acupuncture tip. Consult a doctor if necessary
Stimulation is interrupted during therapy	The acupuncture tip is not positioned correctly	Reposition acupuncture tip
	Contact with the contact electrode has not been established	Ensure touch contact with the contact electrode
	The batteries are empty	Replace with new batteries

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

- If the device is to be used by several persons, the acupuncture tip used must be either sterilised or wiped at least 5 times with isopropyl alcohol or a disinfectant prior to use.
- Do not place the unit in direct contact with water or other fluids.
- Switch the unit off when it is not in use.
- Before putting the **Promed ACP** into storage for a longer period of time, remove the batteries from the compartment. Leaking batteries can damage the unit.
- Keep the unit and its accessories in a cool, dry place in the case provided.
- Handle the electrodes by their edges when removing them. To avoid damaging them, never pull at the electrode leads themselves.
- Do not make any sharp kinks in the connecting leads or electrodes.
- After use, stick the electrode pad onto the protective plastic film.
- Do not expose the device to direct sunlight and protect it against dirt and moisture.
- Never place any heavy objects on the device.
- You can clean your **Promed ACP** by wiping it over carefully with a cloth moistened with mild soapy water. You may also use isopropyl alcohol or a soap solution. Household detergents and cleaning products are not suitable.
- The unit must not be returned to your dealer or the manufacturer for inspection and recalibration. If you need such inspections or recalibrations because of your internal Quality Management System don't hesitate to contact Promed.
- If you should experience other problems, consult your dealer, returning the unit if necessary. Never attempt to repair a fault yourself.

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

**Only for EU countries:** Do not throw power tools in with the domestic waste! In accordance with the European Directive 2012/19/EU on waste electrical and electronic equipment and its implementation into national law, nonfunctioning power tools must be collected separately and recycled in an environmentally suitable manner.

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices using the appropriate collection systems or send the device to the place of purchase for disposal. The place of purchase will then recycle the device.

#### **Rechargeable batteries/batteries:**

Do not throw your rechargeable batteries/batteries in the domestic waste, in fire, or in water. Rechargeable batteries/batteries should be collected, recycled, or disposed of in an environmentally friendly manner.

#### **Only for EU countries:**

According to the directive 91/157/EEC, defect or used rechargeable batteries/batteries must be recycled. Waste rechargeable batteries/batteries can be directly handed to:

**Promed GmbH**, Lindenweg 11, D-82490 Farchant.

This product was produced with the utmost care and was thoroughly tested before it left the factory. Therefore, we provide a warranty period of 24 months from the date of purchase for this product.

**Promed** products conform to their description and the relevant specifications; it is your responsibility to ensure that the product which you have purchased is designed for your intended use.

- In the event of material defects and/or manufacturing errors which appear during proper use and are detected during the warranty period, we will replace all defective parts of the product free of charge within the warranty period, including the labour costs for warranty repairs.
- **The warranty does not cover:**
  - Normal wear on the product
  - Defects which were caused during transportation or storage of the product
  - Defects or damage which were caused through improper use or inadequate maintenance
  - Damage resulting from disregarding the notes in the operating manual
  - Damage resulting from modifications to the product which were not carried out by **Promed**.
  - Damage caused by sharp objects or as a result of torsion, compression, being dropped, an abnormal impact or other actions which are beyond **Promed**'s reasonable control.
  - Parts which are subject to wear (e.g. moving parts such as ball bearings, etc., fasteners, etc.) are generally excluded from the warranty.
- **A claim under the warranty is not enforceable if:**
  - the product is not sent back in its original packaging or in adequately secure packaging;
  - it has been modified or repaired by a person or company other

than **Promed** or one of **Promed**'s authorised dealers;

- the product has been repaired using parts which are not approved by **Promed**;
- the serial number / LOT number of the product has been removed, deleted, changed or made illegible.
- For reasons of hygiene, any products which have been exposed to direct contact with a body or with bodily fluids (e.g. blood) must be packaged in an additional plastic bag before they are sent back. In this case, a specific note of this situation must be made in the package or in the accompanying documents.
- For components or products which are replaced during the warranty period, a warranty is only provided for the remainder of the original warranty period; provided that the product or part replacement was carried out by **Promed** or one of **Promed**'s authorised dealers.
- The warranty period begins on the day of purchase. Claims under the warranty must be established within the warranty period. Complaints arising after the expiry of the warranty period cannot be considered.
  - Under these warranty conditions, the warranty shall only apply if the date of purchase can be proven by means of a sales receipt or similar.
  - We reserve the right to make technical and visual changes as well as changes to the features.
- This warranty is only valid and enforceable in the country in which it was purchased by the original buyer, provided that it was **Promed**'s intention that the product should be offered for sale in this country. This warranty is also enforceable in any country within the European Economic Area in which **Promed** has an authorised importer or sales partner. Depending on each country, specific and differing warranties and guarantees may

DE

EN

FR

IT

ES

NL

PL

RU

TR

AR

be applicable owing to the particular applicable legislation. These legal provisions are neither excluded nor limited by these warranty conditions. To the extent permissible under national law, the warranty period shall not be extended, renewed or otherwise affected by any subsequent resale, repair or replacement of the product.

- The provisions of the UN Convention on Contracts for the International Sale of Goods shall not apply.
- Our warranty conditions do not affect the seller's legal warranty obligation.
- To the fullest extent possible under the applicable mandatory legislation, these warranty conditions constitute your sole and exclusive legal remedy and shall apply in place of any other express or implied warranty conditions. **Promed** is not liable for any unusual, incidental, punitive or consequential damages including, but not limited to loss of profits, loss of use, loss of revenue, costs of replacement equipment or facilities, third party insurance claims, and damage to property which are the result of the purchase or use of the product or occur owing to a breach of warranty, a breach of contract, negligence, product defects or any other legal or statutory circumstances, even though **Promed** was aware of the possibility of such damages. **Promed** is not liable for any delay in warranty claims.
- **Promed** cannot be held liable for any possible translation errors.
- The following information is vital for problem-free processing:
  1. Original sales slip/receipt or dealer's stamp with the date of purchase
  2. Defect detected
  3. Device name / type / serial / lot number

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Transportation only in the original packaging

### Transportation/storage conditions

Transportation and storage temperature:  
-25° C – +70° C  
Relative humidity:  
15 % - 95 %  
non-condensing  
Air pressure:  
700 to 1060 hPa

### Operating conditions:

Operating temperature:  
+5° C – +40° C  
Relative humidity:  
15 % - 93 %  
non-condensing  
Air pressure:  
760 to 1060 hPa

A sudden change in temperature is to be avoided due to the possible condensation of the device. Use the device only after a temperature equalisation.

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